GUÍA RÁPIDA VITEK® 2

VERSION 2020

PIONEERING DIAGNOSTICS

COMPONENTES DEL SISTEMA VITEK® 2 COMPACT

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Work List

Number of Isolates: 24

	Isolate Status	Accession ID #	Organism	ABS Findings	Card Type	Card Status	Test Date	Setup Tech
		1-1	Fluorophora		AST-WT10	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		2-1	V. proteomixia		AST-WT10	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		3-1	Parabac agglutinans		AST-WT10	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		100-1	Staphylococcus		GP	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		104-1	Staphylococcus		GP	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		105-1	Staphylococcus		GP	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		201-1	Staph. aureus		GP	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		302-1	Staph. aureus (MRSA)		AST-0807	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		303-1	Staph. aureus		GP	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		401-1	Proteus mirabilis		GP	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		402-1	Proteus mirabilis		AST-0846	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		502-1	K. pneumoniae		GP	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		503-1	K. pneumoniae		AST-0846	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		504-1	K. pneumoniae		GP	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		606-1	K. oxytoca		AST-0846	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		801-1	Enterobacter		GP	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		901-1	Ent. faecalis		GP	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)
		1002-1	Escherichia		GP	Final	Wed Dec 20 2017	Laboratory Department (Lab Dept)

1 - 20 of 24 Items

User Name: Laboratory Administrator (Lab Admin)

VITEK 2 Systems - Release version: Server: VITEK-0000

- ☒ QC
- ☒ AES Findings
- ☒ Alert
- ☒ Bar Code
- ☒ Card Status
- ☒ Card Type
- ☒ First Name
- ☒ Isolate Status
- ☒ Last Name
- ☒ Lot #
- ☒ McFarland
- ☒ Org Origin
- ☒ Organism
- ☒ Patient ID
- ☒ QC Dev
- ☒ Setup Tech
- ☒ Test Date

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Organism: Staph.saprophyticus

AES Findings: Consistent

Phenotypes Selected for Review: OXAZOLIDINONE RESISTANT

Card Comments:

Advanced Reporting Tool Comments:

AST-GP67!

Antibiotic	MIC	INT	Antibiotic	MIC	INT	Antibiotic	MIC	INT
<input type="checkbox"/> Beta-Lactamase	POS	+	<input type="checkbox"/> Ciprofloxacin	≤0.5	S	<input type="checkbox"/> Vancomycin	1	S
<input type="checkbox"/> Cefoxitin Screen	NEG	-	<input type="checkbox"/> Levofloxacin	0.5	S	<input type="checkbox"/> Tetracycline	≤1	S
<input type="checkbox"/> Benzpenicillin	0.12	R	<input type="checkbox"/> Moxifloxacin	≤0.25	S	<input type="checkbox"/> Tigecycline	≤0.12	S
<input type="checkbox"/> Ampicillin			<input type="checkbox"/> Inducible Clindamycin Resistance	NEG	-	<input type="checkbox"/> Nitrofurantoin	≤16	S
<input type="checkbox"/> Oxacillin	0.5	S	<input type="checkbox"/> Erythromycin	≤0.25	S	<input type="checkbox"/> Rifampicin	≤0.5	S
<input type="checkbox"/> Gentamicin High Level (synergy)			<input type="checkbox"/> Clindamycin	≤0.25	S	<input type="checkbox"/> Trimethoprim/Sulfamethoxazole	≤10	S
<input type="checkbox"/> Streptomycin High Level (synergy)			<input type="checkbox"/> Quinupristin/Dalfopristin	0.5	S			
<input type="checkbox"/> Gentamicin	≤0.5	S	<input checked="" type="checkbox"/> Linezolid	≥8				

Patient Name: Schleicher, Butch

Analysis Status: 14.93 hr - Final

AST Offline Tests: Beta-Lactamase +

Analysis Messages: The following antibiotic(s) are not claimed: Ampicillin, Gentamicin High Level (synergy), Streptomycin High Level (synergy).

McFarland: 0.56

Setup Tech: Laboratory Supervisor (LabSuper)

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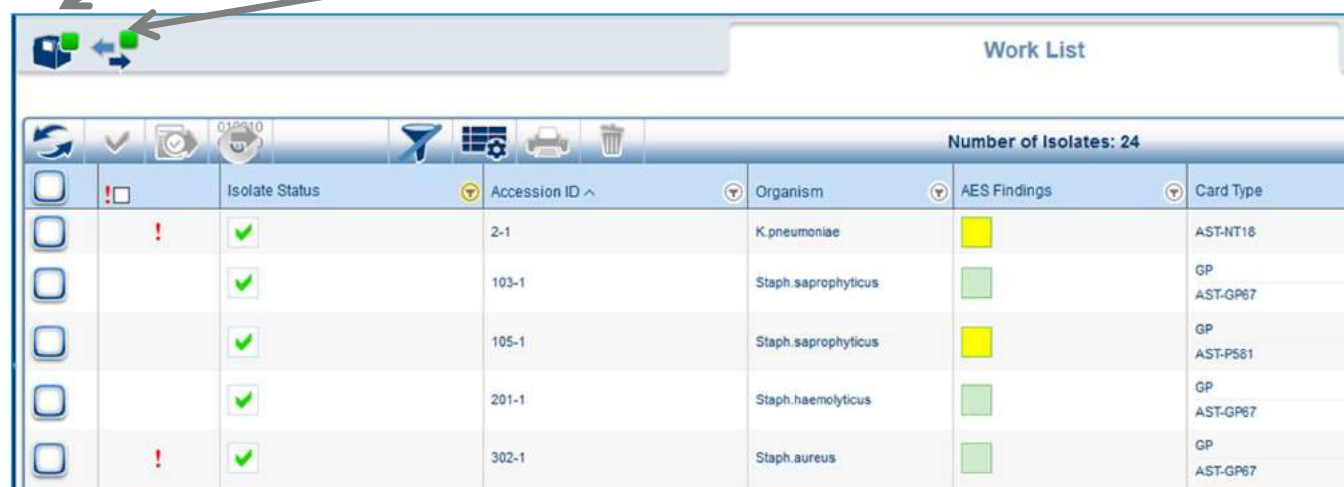
The screenshot displays the VITEK 2 software interface for Quality Control (QC) setup. At the top, there are tabs for 'Sample Preparation' and 'Configuration'. The main window is titled 'VITEK® 2' and features a 'QC' icon circled in green. Below this, there are input fields for 'Lab ID', 'Isolate', 'Card Type' (set to 'AST-GP63'), 'QC Organism ID', 'Shipment Date' (set to 'Apr 28, 2018 19:49 CDT'), 'Lot Number' (set to '1342680293'), and 'Expiration Date' (set to 'Mar 29, 2025 13:00 CDT'). A table on the left lists slots 1 through 10, with slot 1 selected. The bottom of the screen shows the VITEK FLEXPREP Release version 1 and the date 20-01.

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